

# Prophylaxis and Periodontology





# Prophylaxis and Periodontology

Company	4-5
Product overview Piezo Scaler	6-7
Product overview Air scaler, Straight and contra-angle handpieces	8–9
Timon , Diana Caslar	10-11
Tigon+ Piezo Scaler	10-11
Tigon Piezo Scaler	12-13
Pyon 2 Piezo Scaler	14-15
Piezo tips	16-17
Proxeo straight and contra-angle handpieces	18-19
	20 24
Proxeo air scaler	20–21
Technical information	22–23
reciffical information	22-23



W&H Dentalwerk Bürmoos



W&H Sterilization Brusaporto



Healthy development occurs when action is taken in advance. This applies for prophylaxis and periodontology just as it does for W&H: after all, as a family-owned company, we have been operating with the appropriate far-sightedness for generations.

The secret to our success: W&H has always focussed on people — regardless of whether they are partners, customers, patients or employees. It would only be natural for you now to ask critically: Is that really put into practice? What we suggest: Put us to the test!



# People have Priority

#### W&H at a glance

- > Head office in Austria
- > 1000 employees worldwide
- > 80 teams
- > Established in Berlin in 1890 by Jean Weber and Hugo Hampel
- > 3 production facilities
- > 19 subsidiaries in Europe and Asia
- > Export share over 95%

For more information see wh.com



#### Whatever happens:

With W&H prophylactic and periodontic products you are always prepared.



»Today we can say to the patient that tooth scaling is not a painful treatment.«

Dr. Jure Robic, dentist, Slovenia



»Tigon's technology takes everything into account — the result is increased patient comfort and more relaxed working for the dentist.«

Babette Mc Laren-Thomson, dental assistant, Germany



»The delicate characteristics of the tips guarantee a minimally invasive procedure. The Pyon 2 makes easy work of what is often not possible by hand.«

Tania Eberle, dental hygienist, Germany

#### Your advantages at a glance

- > Temperature-controlled liquid for patient well-being
- > Optimally coordinated tip range
- > Handpiece with 5x Ring LED
- > Large, easy-to-read display
- > Self-explanatory menus
- > Multifunctional foot control
- > Tank can be thermo washer disinfected and sterilized
- > Handpiece can be thermo washer disinfected and sterilized

# Tigon+ Piezo Scaler



#### Your advantages at a glance

- > Easy to use
- > Handpiece with 5x Ring LED
- > Coolant tank with 400 ml
- > Optimally coordinated tip range
- > Removable and sterilizable tray
- > Adjustable and removable handpiece holder
- > Perfect electronics
- > Handpiece can be thermo washer disinfected and sterilized

## Tigon Piezo Scaler



#### Your advantages at a glance

- > Handpiece with 5x Ring LED
- > Work with daylight quality
- > Exact reproduction of tooth colour
- > Significantly improved contrast sensitivity
- > Optimally coordinated tip range
- > Handpiece can be thermo washer disinfected and sterilized

### Pyon 2 Piezo Scaler



More information on page 14



»I only use W&H Proxeo products.
They have both a comfortable size and a comfortable surface.«

Lolita Jadermark, dental hygienist, Sweden



»A handy and high-performance sonic scaler with a wide spectrum of uses in periodontology.«

Dr. Christian Graetz, senior dentist at a university hospital, Germany

#### Your advantages at a glance

- > Long life with unique sealing system
- > Best view with small head
- > Efficient treatment with ergonomic design
- Optimal speeds for atraumatic and efficient polishing/cleaning
- > Easy to clean and hygiene-friendly
- > Thermo washer disinfectable and sterilizable



#### Your advantages at a glance

- > Stepless regulation
- > Large optic outlets
- > Efficient range of tips
- > Integrated spray
- > First-class hygiene
- > Thermo washer disinfectable and sterilizable



Proxeo air scaler

More information on page 20

# Tigon+ Piezo Scaler



## The patient's best friend

**Simple, safe and relaxing.** The piezo scaler Tigon+ facilitates your daily work processes considerably and your patient remains completely relaxed during the treatment.



Very gentle thanks to temperaturecontrolled liquid Even sensitive teeth and gums are no longer irritated.



Safe treatment thanks to three preset power modes

Power-Basic-Smooth — you can set the scaler power to suit your preferences.



Five different treatment options
The navigation tooth on the display
shows which of the five programmes
you are currently using.



Handpiece with 5x Ring LED for optimal lighting conditions.



Perfectly coordinated tip range For prophylaxis and periodontology, endodontics, implantology and restoration.



Multifunctional foot control Simple operation so you can concentrate on your patient.

# Tigon Piezo Scaler



## Perfection made easy

#### User-friendly, compact and easy to use.

The Tigon has a range of practical advantages for daily use.



Easy to use

The scaler power and coolant flow rate can be varied using the buttons.



Handpiece with 5x Ring LED

Perfect lighting conditions and the same visual perception as with daylight.



Perfectly coordinated tip range

For prophylaxis and periodontology, endodontics, implantology and restoration.



 $400\,\text{ml}$  coolant tank

Scale with 50 ml intervals so you can mix the rinsing liquids directly in the coolant tank.



Tray and handpiece holder

The tray for your instruments can be thermo washer disinfected and sterilized.



Ergonomically formed handpiece

For safe working without fatigue.

# Pyon 2 Piezo Scaler



## Easy to use, top results

**High-quality, ergonomically designed and safe to use.** The handpiece with five LEDs guarantees fatigue-free working with a considerably larger illuminated field.



Five LEDs
Robust and sterilizable light sources
guarantee optimal illumination.
Significantly longer life than
conventional light sources.



Two versions
The Pyon 2 handpiece is available with a 5x Ring LED or without light.



Large tip range
Best results in endodontics,
prophylaxis, periodontology
and micropreparation.



The ergonomic design of the handpiece and its light weight allows you to work for longer without tiring.



Separable
The Pyon 2 handpiece can be taken apart, is easy to clean, can be sterilized and is thermo washer disinfectable.




**Piezo built-in kit**For installation in your existing dental unit.

# Piezo tips



Tip changer

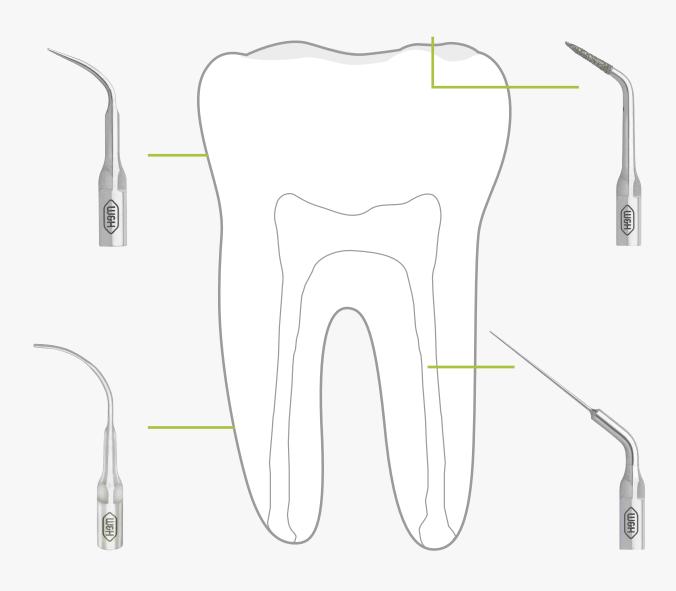



Restoration



Universal



Periodontology



**Endodontics** 

## Compact, user-friendly, safe

Compact, high-quality tip system. 22 universal and special tips, perfectly coordinated for the requirements of the respective application areas and individual treatment steps. All W&H tips can be thermo washer disinfected and sterilized.



Universal tips and implant clean tip 1U, 2U, 3U, 4U, 1I Implant




**Periodontal tips**1P, 3Pr/3Pl, 4P, 5Pr/5Pl, 5P



Restoration tips 1R, 2R, 3Rm, 3Rd



**Endo tips** 1E, 2E, 3E, 4E, 5E, 6E

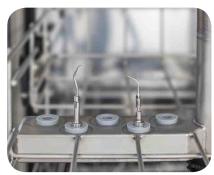

W&H adaptor set
The W&H adaptor set makes it
possible for you to clean and
disinfect the tips automatically
using a validated procedure in the
thermo washer disinfector — inside
and outside.

# Proxeo straight and contra-angle handpieces



## For a beaming smile

**Compact and simple.** Proxeo, the contra-angle handpiece for cleaning, polishing and fluoridation. The most important advantages: its small head with the special sealing system, which prevents any polishing paste from entering.



Long life and hygienic
Long service life with special
sealing system. Prevents
penetration of paste or particles
resulting from treatment.



Best view
A small contra-angle handpiece
head ensures an excellent view
of treatment sites — even in the
distal region.



Incredibly adjustable
The UNIVERSAL contra-angle
handpiece can be used for
all standard screw-in caps,
brushes and Snap-On adaptors.



Versatile and flexible
YOUNG system — can be used for
YOUNG screw-in caps and brushes.
The sealing effect is renewed every
time the cap is changed, which
ensures maximum hygiene.



Simply fast
The WP-66 M contra-angle handpiece
with push-button chucking system
is compatible with all cups,
brushes and Snap-On adaptors
with 2.35 mm shaft.



User-friendly and hygiene-friendly
The straight handpiece for all
disposable contra-angle handpieces
with plastic shaft, Doriot system.
The complete contra-angle
handpiece is removed after use.

## Proxeo air scaler



## Prophylaxis with style

**Fast, gentle, powerful.** The air-driven scaler with the integrated spray and two large optic outlets is ideal for the removal of plaque and calculus as well as cleaning root surfaces. The result: your patients feel better in general and more self-confident thanks to their clean teeth surfaces.



#### Stepless regulation

The powerful, controlled deployment of the power can be adjusted to suit the application in question using the power regulator.



Large optic outlets

Two large optic outlets provide perfect illumination of the treatment site.



Efficient range of tips

The combination of the scaler with the tip range is the ideal basis for perfect results.



Integrated spray

The integrated spray continuously cools the treatment site and simultaneously washes away the abraded microparticles and deposits.



#### First-class hygiene

Their smooth surface means that the instruments are easy to clean and can be both thermo washer disinfected and sterilized.



W&H adaptor set

The W&H adaptor set makes it possible for you to clean and disinfect the tips automatically using a validated procedure in the thermo washer disinfector — inside and outside.

# Tigon+ Piezo Scaler

#### **C €** 0297

Туре:	PA-123 / PA-115	
Mains voltage:	220 - 240 V / 110 - 130 V	
Rated current:	0.1 - 1.6 A / 0.1 - 2.5 A	
Permissible voltage fluctuation:	± 10%	
Frequency:	50 – 60 KHz	
Max. output power to handpiece under load (ultrasonic):	10 W	
Length of handpiece cable:	1.8 m	
Height / width / depth:	178 x 276.5 x 300 mm	
Foot control:	S-N1, C-NF	
Weight:	3.9 kg	

# Tigon Piezo Scaler

#### **C €** 0297

Туре:	PB-323 / PB-315	
Mains voltage:	220 - 240 V / 110 - 130 V	
Rated current:	0.1-0.8A/0.1-1.6A	
Permissible voltage fluctuation:	± 10%	
Frequency:	50 – 60 KHz	
Max. output power to handpiece under load (ultrasonic):	10 W	
Length of handpiece cable:	1.8 m	
Height / width / depth:	178 x 276.5 x 300 mm	
Can be used with foot control:	C-NF	
Weight:	2.7 kg	

# Pyon 2 Piezo Scaler

#### **C €** 0297

Туре:	PB-200 LED	PB-200
Handpiece:	5x Ring LED	without light
Supply hose with coupling:	1.8 m	1.8 m
Handpiece weight:	52.5 g	52.5 g
Coolant hose:	2 m	2 m
Coolant volume:	approx. 50 ml/min approx. 50 ml/min	
Mains voltage:	rage: 100 – 240 V, 50/60 Hz 100 – 240 V, 50/60 Hz	

# Proxeo straight and contra-angle handpieces



Туре:	WP-64 MU	WP-64 M	WP-66 M	HP-44 M	
Reduction ratio:		4:1			
Coupling system:	ISO 3964 (DIN 13.940)				
Applicable attachments:	All standard screw-in caps, brushes and snap-on caps (screw thread 1-72 UNF 2A)	YOUNG Screw-in caps and brushes	Caps, brushes and snap on caps with 2.35 mm shaft	Disposable contra-angle handpiece with plastic shaft, Doriot system	
Max. motor speed:	10,000 min <sup>-1</sup>				
Maximum speed for caps and brushes:	2,500 min <sup>-1</sup>				
Recommended speed for caps and brushes:	1,250 – 1,500 min <sup>-1</sup>				

#### Proxeo air scaler



Type:	ZA-55 L / ZA-55 LM / ZA-55 LS	ZA-55 / ZA-55 M	ZE-55 BC	ZE-55 RM
Light:	yes	no	no	
Power regulation:	yes	no	no	
Oscillating frequency:		6 kHz		
Oscillation amplitude:	< 200 μm			
Operating pressure:	$3 \pm 0.3  \text{bar} /  3^{+0.6}_{-0.3}  \text{bar} /  3 \pm 0.3  \text{bar}$	$3\pm0.3\mathrm{bar}/3^{+0.6}_{-0.3}\mathrm{bar}$	2.2 – 2.5 bar	
Coolant flow volume:	30 – 60 ml/min		min. 30 ml/min	
Coupling:	RQ-24, RQ-34 / RM-34 LED, Multiflex®* / Sirona®*	RQ-3, RQ-4, RQ-14 / Mulitflex®*	lulitflex®* Fixed connection	
Applicable W&H tips:	1AU, 2AU, 3AU, 1AP, 2API, 2APr, 3AP Tips from other manufacturers must be approved for use with W&H air scalers.			
In the scope of delivery:	Tip 1AU, 2AU, 3AU, 1AP, Tip changer		2 x Tip 1AU	
Thermo washer disinfectable:	yes		no	

 $<sup>^{</sup>st}$  Multiflex and Sirona are trademarks of third parties not affiliated with W&H Dentalwerk Bürmoos GmbH.

# Prophylaxis and Periodontology

Scan the QR-Code with your mobile telephon for more information on our products



Pictures are for illustrative purposes only. Additional equipment and accessories shown are not included as standard.

W&H Dentalwerk Bürmoos GmbH Ignaz-Glaser-Straße 53, Postfach 1 5111 Bürmoos, Austria t + 43 6274 6236-0 f + 43 6274 6236-55 office@wh.com wh.com